


Specifications

EndoStratus™ Irrigation Pump
REF#: EGA-500
Dimensions: 200(W) x 173(H) x 385(D)mm
Weight: 11.4 lbs
Power Requirements: 100-240 V
Max Flow Rate:
 Auxiliary water port: 0-280 ml/min
 Biopsy channel: 0-600 ml/min



EndoStratus™ CO₂ Insufflator
REF#: EGA-501
Dimensions: 200(W) x 173(H) x 385(D)mm
Weight: 12.5 lbs
Power Requirements: 100-240 V

Easy Set Up and Installation

1 Choose a CO₂ Source

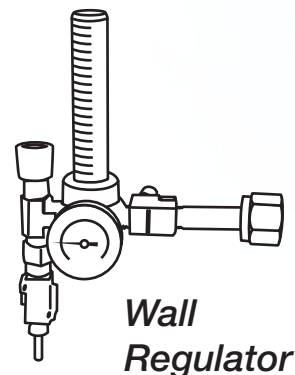

2 Connect Hose and Tubing to Insufflator

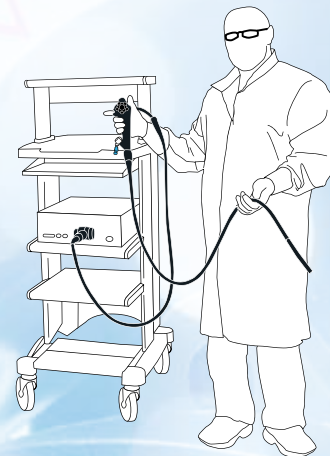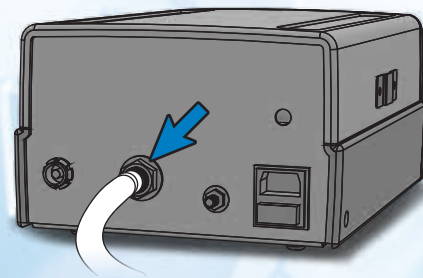

3 Ready for Procedure

Compatible Accessories

Tubing



REF #100130
EndoGator™ Tubing

REF # 100145CO2
Endo SmartCap™ Tubing for Olympus® 140/160/180 Series Endoscopes

REF # 100165CO2
Endo SmartCap™ Tubing for Fujinon® Endoscopes and Pentax® Endoscopes

REF # 100150CO2
Endo SmartCap™ Tubing for Olympus® 100/130 Series Endoscopes and Pentax® Endoscopes



REF #100603
EndoGator® Kit Tubing

REF #100611
EndoGator® Kit Tubing with CO₂

REF #100609
EndoGator® Hybrid Tubing with CO₂

Compatible with Olympus® Endoscopes

Connectors



REF #100241
Single Use Auxiliary Water Jet Connector for Olympus® Endoscopes

REF #100242
Single Use Auxiliary Water Jet Connector for Pentax® Endoscopes

REF #100141
Endo SmartCap® Adapter for Fujinon® 530 Series Colonoscopes

Biopsy Irrigation



REF #100303
DEFENDO® Y-Opsy® Irrigator
for Olympus® or Fujinon® Endoscopes

MEDIVATORS®
ENDO stratus™
 IRRIGATION PUMP and CO₂ INSUFFLATOR



MEDIVATORS®, EndoStratus™, EndoGator®, Endo SmartCap®, DEFENDO®, and Y-Opsy® are trademarks or registered trademarks of Medivators Inc.
 Fujinon® is a registered trademark of Fuji Photo Film Co., Ltd. and/or its affiliated entities.
 Olympus® is a registered trademark of Olympus Corporation, Olympus America, Inc., and/or its affiliated entities.
 PENTAX® is a registered trademark of PENTAX of America, Inc.

Distribuidor

Beracah MEDICA

ventas@beracahmedica.mx

662 15 10 979

MEDIVATORS®
 A Cantel Medical Company

The CO₂ Solution for GI Endoscopy \\\\\\\

The MEDIVATORS® **EndoStratus™** CO₂ Insufflator is a versatile system that connects to either a wall or tank CO₂ source. Combined with the full line of the disposable MEDIVATORS® **Endo SmartCap®** CO₂ Tubing, the **EndoStratus™** CO₂ Insufflator offered by MEDIVATORS® allows clinicians to offer the optimal procedure for the patient.

MEDIVATORS® offers a complete range of CO₂ solutions to suit individual facility needs for *Olympus®*, *Pentax®* and *Fujinon®* GI endoscopes. With direct connection to sterile water and any CO₂ tank or wall source - MEDIVATORS® has the complete solution for your facility and patients.

- CO₂ is easier on the patient both during and after the procedure due to rapid absorption by surrounding tissue.¹
- Minimizes patient discomfort associated with bloating and cramping from room air.¹



Daily Use *Endo SmartCap*® CO₂ Source Tubing

Warm Water Technology


Deliver irrigation for GI endoscopy procedures with a warm water option by utilizing the **EndoStratus™** Irrigation Pump. Combine the pump with any GI endoscope and irrigate through the auxiliary water jet or biopsy channel to improve procedural visibility.

- Greater visibility for improved diagnosis and polyp detection
 - Auxiliary water jet or biopsy channel irrigation capability
 - Warm water technology
 - Compatible with MEDIVATORS®
- EndoGator®** Tubing and Connectors
- Easy to use water bottle cradle
 - Automatic prime button provides instant irrigation upon foot pedal depression



Auxiliary Water Jet Connection


CO₂ Insufflator Features

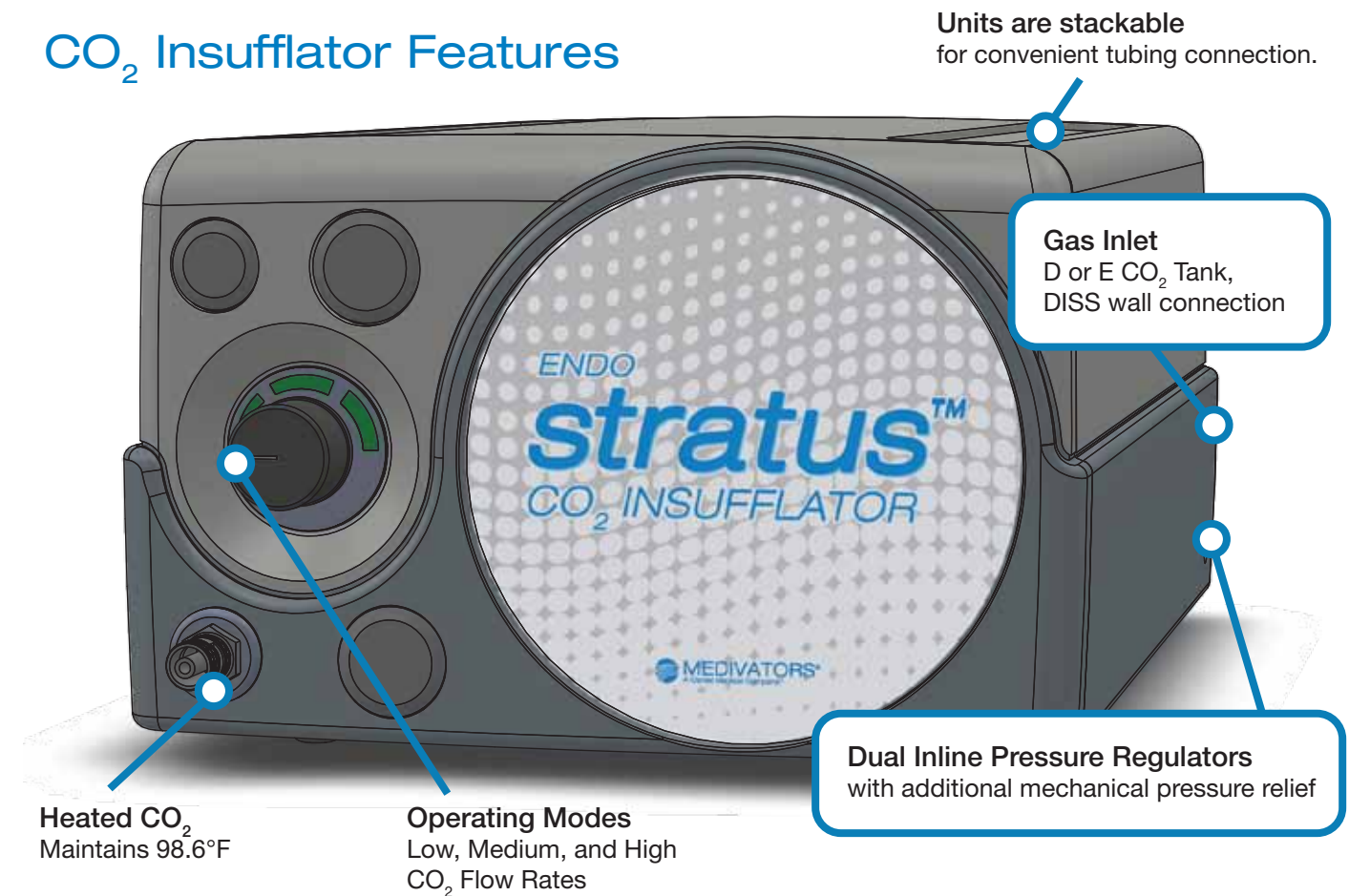

Irrigation Pump Features

Safety Start Feature

The unit will not operate unless the pump head is closed.

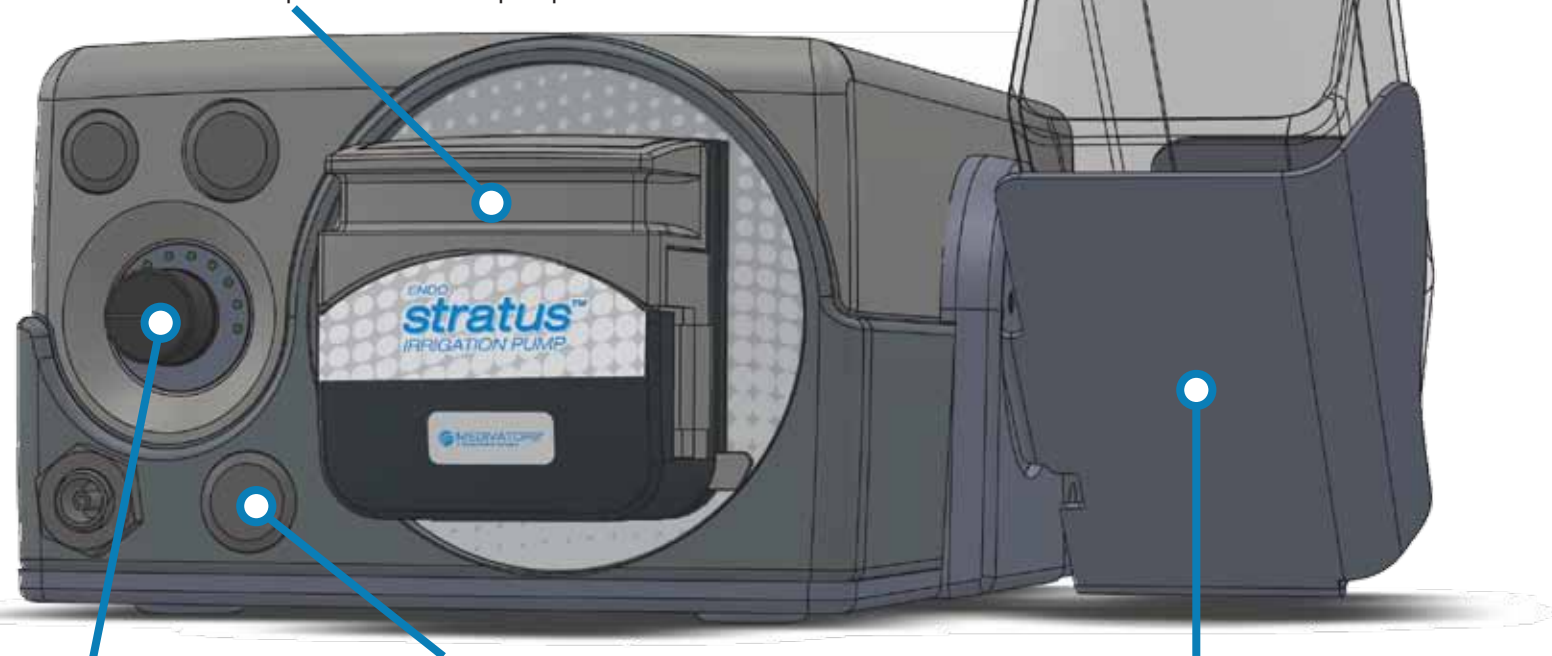

Flow Rate
Auxiliary water port: 0-280 ml/min
Biopsy channel: 0-600 ml/min

Prime Button
Runs the pump for 20 seconds allowing the user to prime the tube set without pressing foot pedal.

Water Warming Angled Bottle Holder
Maintains 98.6°F water temperature.